

RADDOSE-3D Command Reference

Dated: 26 July 2024

Please cite the following primary publication when using RADDOSE-3D:

RADDOSE-3D: time- and space-resolved modelling of dose in macromolecular crystallography

Zeldin, O. B.; Gerstel, M. & Garman, E. F. (2013). *J. Appl. Cryst.* **46**
doi:10.1107/S0021889813011461

For intensity decay models or doses for electron diffraction (RADDOSE-ED), please cite:

Doses for X-ray and electron diffraction: New features in RADDOSE-3D including intensity decay models

Dickerson J. L.; McCubbin P. T. N.; Brooks-Bartlett, J. C., & Garman, E. F. (2024). *Protein Science* **33**
doi.org/10.1002/pro.5005

If calculating doses for X-ray Free-Electron Lasers (RADDOSE-XFEL), please cite:

RADDOSE-XFEL: femtosecond time-resolved dose estimates for macromolecular X-ray free-electron laser experiments

Dickerson J. L.; McCubbin P. T. N., & Garman, E. F. (2020). *J. Appl. Cryst.* **53**
doi.org/10.1107/S1600576720000643

If using microcrystals or microbeams, we recommend using the Monte Carlo simulations outlined in:

Doses for experiments with microbeams and microcrystals: Monte Carlo simulations in RADDOSE-3D

Dickerson J. L. & Garman, E. F. (2020). *Protein Science* **30**
doi.org/10.1002/pro.3922

For dose calculations related to SAXS experiments, please cite:

Development of tools to automate quantitative analysis of radiation damage in SAXS experiments

Brooks-Bartlett, J. C.; Batters, R. A.; Bury, C. S.; Lowe, E. D.; Ginn, H. M.; Round, A., & Garman, E. F. (2017). *J. Synch. Rad.* **24**
doi:10.1107/S1600577516015083

For developments in RADDOSE-3D, including the ability to incorporate experimental profiles and irregular polyhedral crystal geometries, as well as the inclusion of models for energy loss from the crystal volume due to the Compton effect, please additionally refer to:

Estimate your dose: RADDOSE-3D

Bury, C. S.; Brooks-Bartlett, J. C.; Walsh, S. P. & Garman, E. F. (2018). *Protein Science* **27**
doi:10.1002/pro.3302

If there are heavy atoms present in the sample, it may be more accurate to simulate fluorescent escape, with the implementation detailed in:

Structural and Electronic Effects of X-ray Irradiation on Prototypical [M(COD)Cl]₂ Catalysts

Fernando, N. K. et al. (2021). *The Journal of Physical Chemistry A* **125**

doi:10.1021/acs.jpca.1c05759

Contents

1	General syntax considerations	4
2	Crystal block	4
2.1	TYPE	4
2.2	WIREFRAMETYPE	5
2.3	MODELFILE	5
2.4	DIMENSION	5
2.5	PIXELSPERMICRON	6
2.6	ANGLEP	6
2.7	ANGLEL	6
2.8	CONTAINERMATERIALTYPE	7
2.9	MATERIALMIXTURE	7
2.10	MATERIALELEMENTS	7
2.11	CONTAINERTHICKNESS	7
2.12	CONTAINERDENSITY	8
2.13	ABSCOEFCALC	9
2.14	PDB	10
2.15	SEQFILE	10
2.16	CIF	11
2.17	UNITCELL	11
2.18	NUMMONOMERS	11
2.19	NUMRESIDUES	11
2.20	NUMRNA	11
2.21	NUMDNA	12
2.22	NUMCARB	12
2.23	PROTEINHEAVYATOMS	12
2.24	SOLVENTHEAVYCONC	12
2.25	SOLVENTFRACTION	13
2.26	PROTEINCONC	13
2.27	SMALLMOLEATOMS	13
2.28	CALCULATEPEESCAPE	13
2.29	CALCULATEFLESCAPE	13
2.30	CALCSURROUNDING	13
2.31	SURROUNDINGHEAVYCONC	14
2.32	GONIOMETERAXIS	14
2.33	POLARISATIONDIRECTION	14
2.34	DENSITYBASED	14
2.35	SURROUNDINGELEMENTS	14
2.36	SURROUNDINGDENSITY	14
2.37	DDM	14
2.38	SUBPROGRAM	15

3	Beam block	16
3.1	TYPE	16
3.2	FLUX	16
3.3	FWHM	16
3.4	ENERGY	16
3.5	ENERGYFWHM	16
3.6	FILE	17
3.7	PIXELSIZE	17
3.8	COLLIMATION	17
3.9	PULSEENERGY	17
4	Wedge block	17
4.1	EXPOSURETIME	17
4.2	ANGULARRESOLUTION	18
4.3	STARTOFFSET	18
4.4	TRANSLATEPERDEGREE	18
4.5	ROTAXBEAMOFFSET	18

RADDOSE-3D can take input from one or more files and/or from standard input (STDIN). Any input will be processed by the `InputParser` class and the RADDOSE-3D ANTLR parser. This section describes the syntax of accepted input. Advanced users of RADDOSE-3D can create their own input method that need not rely on the `InputParser` class or the RADDOSE-3D ANTLR parser. This feature will not be covered in this reference.

The simplest use case of RADDOSE-3D will involve only one file describing the entire experiment. In some instances it may be desired to split up the input into a number of files, e.g. one file describing the crystal, one automatically updated file describing the current beam on the beamline, and one file chosen from a set of possible wedge strategies. Each file can contain an arbitrary number (including none) of `Crystal`, `Beam` and `Wedge` block (henceforth called *blocks*). However, splitting up blocks across multiple files is not allowed.

The parser will read the input sequentially, and, when multiple sources are given, one source after the other in the specified order. While the parser may accept `Crystal`, `Beam` and `Wedge` blocks in any order, the exposure of a wedge can only take place if both the crystal and the beam have been set either in an earlier file or before the `Wedge` block within the same file.

1 General syntax considerations

Any keywords specified below are case-insensitive. Upper (**CRYSTAL**), lower (**crystal**) and mixed case (**CrYsTaL**) are equivalent.

The characters `#`, `!` and the character sequence `//` denote the start of a comment. Any text from that position until the end of the current line is ignored.

Tabular and newline characters are treated as white space. They can therefore be freely used to format the file for increased readability.

The order of statements within a `Crystal`, `Beam` and `Wedge` block generally is not relevant. There are two exceptions to this rule: The leading keyword (**CRYSTAL**, **BEAM**, **WEDGE**) must be the first keyword of the block. If a keyword is repeated within the same block, then the latter will always override the former.

Every block must be self-contained, e.g. the energy set for the previous `Beam` is not remembered when setting up the following `Beam`, and must be repeated.

Numeric values can be given in scientific notation ($2.0e2 = 2e+2 = 200$), negative values may not have a space between the sign (`'-'`) and the value ($-1.9e-1 = -.19 = -0.19$).

2 Crystal block

A `Crystal` block must begin with the keyword **CRYSTAL**. At least the **TYPE** and **DIMENSION** must be specified. Depending on the chosen **TYPE** further declarations may be required.

2.1 TYPE

With the keyword **TYPE** the underlying crystal implementation is chosen. Currently four distinct crystal implementations exist:

TYPE CUBOID defines a solid crystal with a cuboid shape.

TYPE SPHERICAL defines a solid crystal with a spherical shape.

TYPE CYLINDER defines a solid crystal with a cylindrical shape.

TYPE POLYHEDRON defines an arbitrary crystal shape as a polyhedron. The wire frame type (**WIREFRAMETYPE** keyword - section section 2.2) and file containing the model (**MODELFILE** keyword - section section 2.3) must be specified.

2.2 WIREFRAMETYPE

WIREFRAMETYPE Specifies the type of wire frame model used to model a crystal of polyhedron type. Currently only .obj (geometry definition) files can be used:

WIREFRAMETYPE OBJ specifies the wire frame model to be in .obj format.

2.3 MODELFILE

MODELFILE Specifies the location of the file that contains the wire frame model of the polyhedron crystal. Currently only .obj (geometry definition) files can be read. The models and the .obj files can be generated using the free and open source 3D animation software **BLENDER**. (NOTE: if you are exporting a *Wavefront (.obj)* file in **BLENDER**, then select the option *Triangulate Faces* before you finalise the export. RADDOSÉ-3D works with triangular faces for the polygons.):

MODELFILE SOMEMODELFILE.OBJ specifies the location of the .obj file to be imported.

2.4 DIMENSION

DIMENSION specifies the size of the crystal. Dimensions are given in micrometres (μm). The keyword **DIMENSION** can take either one, two or three parameters:

DIMENSION D with a single number (see section 1) as the parameter *D* is used for specifying the crystal dimensions for spherical crystals (**TYPE SPHERICAL**). The parameter sets the crystal diameter. This syntax cannot be used for cuboid or cylindrical crystals.

DIMENSION D H with two numbers as parameters *D* and *H* is used for specifying the crystal dimensions for cylindrical crystals (**TYPE CYLINDER**). The parameter *D* sets the diameter of the circular cross section of the cylinder. The parameter *H* sets the height of the cylinder. This syntax cannot be used for cuboid crystals.

DIMENSION X Y Z with three numbers as parameters *X*, *Y* and *Z* is used to set the dimensions for cuboid crystals (**TYPE CUBOID**). *X* defines the length of the crystal orthogonal to both the beam and the goniometer at $L=P=0$, (see below) *Y* defines the length along the goniometer axis at $L=P=0$ and *Z* defines the length along the beam axis.

If three parameters are given for a spherical crystal (**TYPE SPHERICAL**) the value for X sets the diameter of the crystal while the values of Y and Z are ignored. If three parameters are given for a cylindrical crystal (**TYPE CYLINDER**) the value for X sets the diameter of the circular cross section, the value for Y sets the height of the cylindrical crystal while the value Z is ignored.

2.5 PIXELSPERMICRON

PIXELSPERMICRON F specifies the resolution of the voxel grid used to represent the crystal in voxels/ μm . Defaults to 0.5 voxels/ μm .

Note: When running RADDOS-3D with large dimension values (with SAXS the capillary dimensions typically range in the millimetre range) RADDOS-3D can terminate with the following error:

Error during invocation of se.raddo.raddose3D.CrystalCylinder: Java heap space. This is because there is too much memory being used. reducing the pixels per micron will sort this problem and also improve the speed of the program. The cost for the speed improvement is the resolution of the sample voxelization.

2.6 ANGLEP

ANGLEP F sets the angle in the plane of the loop between the crystal Y axis and the goniometer axis. The angle is to be given in degrees, but without the degree symbol ($^\circ$). The default P ('plane') angle is 0° .

The rotation angle to be applied to the crystal in the plane of the loop (right handed rotation about Z axis applied to all voxels, as shown in figure 1).

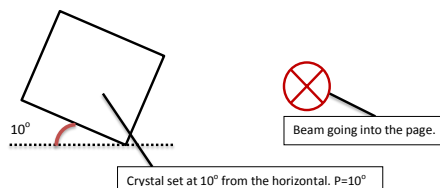


Figure 1: Schematic of **ANGLEP**. Figure courtesy of John Bremridge.

2.7 ANGLEL

ANGLEL F sets the loop angle between the plane of the crystal loop and the goniometer axis. The angle is to be given in degrees, but without the degree symbol ($^\circ$). The default L ('loop') angle is 0° .

The rotation angle to be applied to the angle of the crystal in the loop (right handed rotation about X axis applied to all voxels, as shown in figure 2).

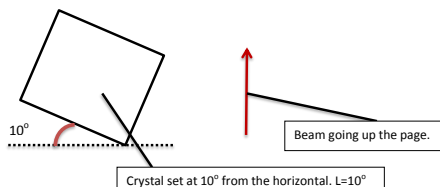


Figure 2: Schematic of **ANGLEL**. Figure courtesy of John Bremridge.

2.8 CONTAINERMATERIALTYPE

CONTAINERMATERIALTYPE MATERIALTYPE Specifies the material type of the container that encases the irradiated sample. This keyword should be used when the sample is encased within a container e.g. SAXS experiment where the sample is irradiated through a capillary. Currently three distinct container material type implementations exist:

CONTAINERMATERIALTYPE NONE The sample is not encased inside any type of container hence there is no attenuation of the beam prior to making contact with the sample. This is the case with a standard X-ray crystallography experiment. This option is selected as the default if the **CONTAINERMATERIALTYPE** is not defined.

CONTAINERMATERIALTYPE MIXTURE Defines a container encasing the irradiated sample which is a mixture of elements, determined by the name of the mixture. If this option is used then the user must specify the material mixture (see section **MATERIALMIXTURE**).

CONTAINERMATERIALTYPE ELEMENTAL Defines a container encasing the irradiated sample in terms of its component elements. If this option is used then the user must specify the list of the material's component elements (see section **MATERIAELEMENTS**).

2.9 MATERIALMIXTURE

MATERIALMIXTURE MATERIAL Specifies the material of the container that encases the irradiated sample by its mixture name. The exact input for **MATERIAL** is given in the URL of the corresponding material on the National Institute of Standards and Technology (NIST) Table 4 webpage (click here). For example, if the material through which the sample is irradiated is pyrex glass then click on the 'Glass, Borosilicate "Pyrex"' link and the URL is:

<http://physics.nist.gov/PhysRefData/XrayMassCoef/ComTab/pyrex.html>.

The material name is given immediately before the ".html" and after the last forward slash "/", in this case the word is "pyrex". Hence to set the container material to pyrex glass the input would be **MATERIALMIXTURE pyrex**. Note that the input is case sensitive and an internet connection is required for this option to work.

2.10 MATERIAELEMENTS

MATERIAELEMENTS E1 I (E1 I (E1 I ...)) Specifies the material of the container that encases the irradiated sample by a list of the component elements of the material. The mass attenuation coefficients for each element are downloaded from the National Institute of Standards and Technology (NIST) Table 3 webpage (click here). For example, if the material through which the sample is irradiated is quartz, which has the formula SiO₂ then the input would be **MATERIALMIXTURE Si 1 O 2**. Note that an internet connection is required for this option to work.

2.11 CONTAINERTHICKNESS

CONTAINERTHICKNESS F Specifies the thickness of the container encasing the irradiated sample. The thickness should be given in μm .

2.12 CONTAINERDENSITY

CONTAINERDENSITY *F* Specifies the density of the container encasing the irradiated sample. The density should be given in grams/centimetre³.

2.13 ABSCOEFCALC

This keyword specifies whether the program should use average absorption and attenuation coefficients, or whether it should calculate them from input crystal parameters.

ABSCOEFCALC AVERAGE

ABSCOEFCALC DUMMY

These two commands are equivalent. Each will cause RADDOS-3D to assume an absorption coefficient of 0.237 mm^{-1} and an attenuation coefficient of 0.281 mm^{-1} . These values are representative of an average crystal at an incident X-ray beam energy of 12.4 keV (1Å). Please see Section 3 in the main paper for more details. Crystal composition keywords will have no effect.

ABSCOEFCALC RD

ABSCOEFCALC RDV2

ABSCOEFCALC RDV3

These three commands are equivalent. RADDOS-3D will call a previous version of RADDOS to estimate absorption and attenuation coefficients.

ABSCOEFCALC RD3D

This command will use the current RADDOS-3D code to calculate the absorption and attenuation coefficients using the crystal composition specified by the user. (NOTE: If the crystal composition is specified then the **RD3D** keyword is preferred over the other **ABSCOEFCALC** keywords, particularly if the crystal contains DNA or RNA molecules).

The composition of the crystal has to be described using the keywords **UNITCELL**, **NUMMONOMERS**, **NUMRESIDUES**, **NUMRNA**, **NUMDNA**, **PROTEINHEAVYATOMS**, **SOLVENTHEAVYCONC** and **SOLVENTFRACTION**. The use of these keywords is described in the sections 2.17–2.38 below. Note that the **SOLVENTFRACTION** keyword is now optional. If it's not given then RADDOS-3D will calculate it.

ABSCOEFCALC EXP

This command should be used if the crystal composition from a Protein Data Bank (PDB) entry is to be used instead of a user specified crystal composition. The current RADDOS-3D code is then used to calculate the absorption and attenuation coefficients using the crystal composition from the PDB entry. The PDB entry must be specified using the **PDB** keyword (section 2.14).

ABSCOEFCALC SEQUENCE

This command should be used if the sample composition from a sequence file is to be used instead of a user specified crystal composition. The sequence file should be in **FASTA** file format. If this command is used then the composition of the crystal can be specified using the following keywords **UNITCELL**, **SEQFILE**, **NUMMONOMERS**, **SOLVENTHEAVYCONC**

ABSCOEFCALC SAXS

This command will use the current RADDOS-3D code to calculate the absorption and attenuation coefficients using the sample composition specified by the user similarly to the keyword **RD3D**. The major difference is that instead of supplying the number of monomers using the keyword **NUMMONOMERS**, the user supplies the protein concentration used in the SAXS experiment using the keyword **PROTEINCONC** (see section **PROTEINCONC**).

Note that RADDOS-3D treats the irradiated samples as stationary objects.

Therefore if the SAXS solution is flowed through the capillary during the exposure then a suitable exposure time for the relevant volume will need to be taken into account.

ABSCOEFCALC SAXSSEQ

This command will use the current RADDPOSE-3D code to calculate the absorption and attenuation coefficients using the sample composition specified by the user similarly to the keyword **SAXS**. The major difference is that you no longer need to give the number of residues or the protein heavy atom composition. Instead the user needs to specify where the sequence file is located using the **SEQFILE**.

ABSCOEFCALC SMALLMOLE

This command will use the current RADDPOSE-3D code to calculate the absorption and attenuation coefficients using the crystal composition specified by the user, similarly to the keyword **RD3D** above. However, this indicates that a small molecule experiment is present, and as such the user must instead specify the composition of the crystal using the keywords **SMALLMOLEATOMS**, **UNITCELL** and **NUMMONOMERS** to explicitly define the entire contents of the small molecule unit cell (including solvent atoms). Note that in this current version, no additional solvent atoms can be specified with **SOLVENTHEAVYCONC** keyword provided below.

ABSCOEFCALC CIF

This option is also provided for small molecule experiments. This command should be used if the crystal composition from a supplied CIF-format file is to be used instead of a user specified crystal composition. The current RADDPOSE-3D code is then used to calculate the absorption and attenuation coefficients using the crystal composition from the CIF entry. A full file path to the CIF-format file must be specified using the additional **CIF** keyword (section 2.16).

2.14 PDB

PDB CODE

Where **CODE** is the four letter code of the PDB entry to be downloaded. From the PDB entry RADDPOSE-3D will calculate **UNITCELL**, **NUMMONOMERS**, **NUMRESIDUES**, **NUMRNA**, **NUMDNA**, **PROTEINHEAVYATOMS** and **SOLVENTFRACTION**, and hence these keywords are not needed. It will NOT calculate **SOLVENTHEAVYCONC** so if there are any heavy atoms in the solvent then they will have to be specified. Note that an internet connection is required for this option to work.

PDB FILENAME

Where **FILENAME** specifies the location of the PDB file for the chosen sample. Only sequence files in the **PDB** format can be read.

2.15 SEQFILE

SEQFILE FILENAME specifies the location of the file that contains the sequence of the chosen sample. Only sequence files in the **FASTA** format can be read.

2.16 CIF

CIF *FILENAME* specifies the location of the file that contains the sequence of the chosen sample. This is exclusively for small molecule experiments and only files in the **CIF** format can be read.

2.17 UNITCELL

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDPOSE (see section 2.13).

UNITCELL A B C

UNITCELL A B C α β γ

Dimensions and angles of the unit cell a , b , c , α , β , γ

The first three numbers specify the unit cell size in Angstroms. The second three numbers optionally specify the unit cell angles alpha, beta and gamma.

The (optional) angles are to be given in degrees, but without the degree symbol ($^{\circ}$). If no angles are specified RADDPOSE-3D assumes default angles of 90° .

2.18 NUMMONOMERS

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDPOSE (see section 2.13).

NUMMONOMERS I specifies the number of monomers in the unit cell. Only integer numbers I should be used. This number should not be confused with the number of monomers in the asymmetric unit.

2.19 NUMRESIDUES

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDPOSE (see section 2.13).

NUMRESIDUES I specifies the number of amino acid residues per monomer. Only integer numbers I should be used. Using this keyword the number and types of atoms are calculated according to the formula

$$\text{amino acid} = 5C + 1.35N + 1.5O + 8H$$

Sulfur atoms, e.g. from CYS and MET residues, should be added explicitly with the **PROTEINHEAVYATOMS** keyword.

The default value for I is 0.

2.20 NUMRNA

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDPOSE (see section 2.13).

NUMRNA I specifies the number of RNA nucleotides per monomer. Only integer numbers I should be used. Using this keyword the number and types of atoms are calculated assuming an average nucleotide content defined as

$$\text{mean nucleotide} = 11.25H + 9.5C + 3.75N + 7O + 1P$$

If a more accurate estimate is required, individual atoms may be entered explicitly with the **PROTEINHEAVYATOMS** keyword.

The default value for *I* is 0.

2.21 NUMDNA

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDose (see section 2.13).

NUMDNA *I*

specifies the number of DNA deoxynucleotides per monomer. Only integer numbers *I* should be used. Using this keyword the number and types of atoms are calculated assuming an average deoxynucleotide content defined as

$$\text{mean nucleotide} = 11.75H + 9.75C + 4N + 6O + 1P$$

If a more accurate estimate is required, individual atoms may be entered explicitly with the **PROTEINHEAVYATOMS** keyword.

The default value for *I* is 0.

2.22 NUMCARB

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDose (see section 2.13).

NUMCARB *I*

specifies the number of carbohydrate residues per monomer. Only integer numbers *I* should be used. Using this keyword the number and types of atoms are calculated assuming an average carbohydrate content defined as

$$\text{mean sugar} = 11H + 6C + 5O$$

If a more accurate estimate is required, individual atoms may be entered explicitly with the **PROTEINHEAVYATOMS** keyword.

The default value for *I* is 0.

2.23 PROTEINHEAVYATOMS

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDose (see section 2.13).

PROTEINHEAVYATOMS *E1 I (E1 I (E1 I ...))* defines a list of atoms to add to the protein part of the absorption. Each species is defined by a two character string *E1* for the elemental symbol, and an integer number *I* of atoms of that species per monomer.

The command **PROTEINHEAVYATOMS S 10 Se 2** would add 10 sulfur and 2 selenium atoms per monomer.

2.24 SOLVENTHEAVYCONC

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDose (see section 2.13).

SOLVENTHEAVYCONC *E1 I (E1 I (E1 I ...))* defines the concentration of elements (not including water) in the solvent in millimoles per litre. Oxygen and lighter elements should not be specified.

The command **SOLVENTHEAVYCONC Na 1000 Cl 1000** specifies 1M sodium chloride in the solvent.

2.25 SOLVENTFRACTION

This keyword only has an effect when the absorption and attenuation coefficients are estimated using a legacy version of RADDOS (see section 2.13).

SOLVENTFRACTION *F*

The fraction of the unit cell that is occupied by solvent. If not given explicitly, this value is estimated from **NUMRESIDUES**, **NUMRNA** and **NUMDNA** using 1.35 g/ml for protein, and 1.35 g/ml for DNA and 1.3 g/ml for RNA.

2.26 PROTEINCONC

This keyword only has an effect when the absorption and attenuation coefficients are estimated using the keyword **SAXS** section 2.13

PROTEINCONC *F*

The protein concentration used in the SAXS experiment in grams per litre.

2.27 SMALLMOLEATOMS

This keyword only has an effect when the absorption and attenuation coefficients are estimated using the keyword **SMALLMOLE** section 2.13

and is only to be used for small molecule experiments

SMALLMOLEATOMS E1 I (E1 I (E1 I . .)) defines a list of atoms to explicitly define the full contents of the unit cell. Each species is defined by a two character string **E1** for the elemental symbol, and an integer number **I** of atoms of that species per ‘monomer’ in the unit cell (note that if the **NUMMONOMERS** keyword is omitted from the input file, then 1 monomer is assumed per unit cell).

The command **SMALLMOLEATOMS C 18 H 15 Bi 8** would add 18 carbon 15 hydrogen and 8 bismuth atoms per ‘monomer’ in the unit cell.

2.28 CALCULATEPEESCAPE

CALCULATEPEESCAPE BOOLEAN Whether or not to include an explicit 3D model for photoelectron escape from the crystal volume. Set to **TRUE** to include or **FALSE** otherwise (defaults to **FALSE** if not set).

2.29 CALCULATEFLESCAPE

CALCULATEFLESCAPE BOOLEAN Whether or not to include an explicit 3D model for fluorescent X-ray escape from the crystal volume. Set to **TRUE** to include or **FALSE** otherwise (defaults to **FALSE** if not set).

2.30 CALCSURROUNDING

CALCSURROUNDING BOOLEAN Whether or not to include photoelectron entry from the surrounding material into the crystal volume. Set to **TRUE** to include or **FALSE** otherwise (defaults to **FALSE** if not set). Note that **CALCULATEPEESCAPE** must be set to **TRUE** for this option to be used.

2.31 SURROUNDINGHEAVYCONC

SURROUNDINGHEAVYCONC *E1 I (E1 I (E1 I ...))* defines the concentration of elements (not including water) in the solvent in millimoles per litre. Oxygen and lighter elements should not be specified.

The command **SURROUNDINGHEAVYCONC Na 1000 Cl 1000** specifies 1M sodium chloride in the surrounding.

2.32 GONIOMETERAXIS

GONIOMETERAXIS *I* specifies the orientation of the goniometer. Currently only 0 and 90 are accepted values, with 0 corresponding to a horizontal goniometer (parallel to the floor) and 90 a vertical goniometer (perpendicular to the floor).

The default value for *I* is 0 (horizontal).

2.33 POLARISATIONDIRECTION

POLARISATIONDIRECTION *I* specifies the direction of the beam polarisation vector. Currently only 0 and 90 are accepted values, with 0 corresponding to a horizontally polarised beam (parallel to the floor) and 90 a vertically polarised beam (perpendicular to the floor).

The default value for *I* is 0 (horizontally polarised).

2.34 DENSITYBASED

DENSITYBASED *BOOLEAN* Whether or not the surrounding is to be defined by elemental composition and density. Set to **TRUE** to indicate density based or **FALSE** otherwise (defaults to **FALSE** if not set).

2.35 SURROUNDINGELEMENTS

SURROUNDINGELEMENTS *E1 I (E1 I (E1 I ...))* defines the formula of a density defined surrounding material. Only to be used when **DENSITYBASED** and **CALCSURROUNDING** are set to **TRUE**

The command **SURROUNDINGELEMENTS C 3 H 8** specifies elemental composition of 3 carbon atoms for every 8 hydrogen atoms

2.36 SURROUNDINGDENSITY

SURROUNDINGDENSITY *F* defines the density of a density defined surrounding material in g/ml. Only to be used when **DENSITYBASED** and **CALCSURROUNDING** are set to **TRUE**

2.37 DDM

DDM *S* specifies the dose decay model used by RADDose-3D. If no DDM is specified, no intensity decay is modelled by default.

DDM LEAL uses the Leal model (Leal RM, Bourenkov G, Russi S, Popov AN. J Synchrotron Radiat. 2013 Jan;20(Pt 1):14-22) which requires three parameters: **DECAYPARAM F F F** which are the γ , B_0 , and β parameters, respectively.

2.38 SUBPROGRAM

This keyword specifies whether the program should use Monte Carlo simulations or RADDOSE-XFEL.

SUBPROGRAM XFEL

Specifies that RADDOSE-XFEL, which calculates the time-resolved dose absorbed by a crystal in an XFEL experiment, should be used.

The pulse energy of the beam must be defined in the beam block with the **PULSEENERGY** keyword (section 3.9).

The pulse length is defined in the wedge block with the **EXPOSURETIME** keyword (section 4.1).

SUBPROGRAM MONTECARLO

Specifies that the RADDOSE-3D Monte Carlo simulations should be run. The thickness of the surrounding material can also be specified in the Monte Carlo simulations using the input **SURROUNDINGTHICKNESS X Y Z**

Both of these subprograms require at least two further inputs:

RUNS *F*

Specifies the number of times the simulation is to be run

SIMPHOTONS *F*

Specifies the number of photons to use in the simulations.

SUBPROGRAM EMED

Specifies that the RADDOSE-ED subprogram should be run, which calculates doses for MicroED experiments.

In the beam block, a further input is required:

EXPOSURE *F*

Specifies the total fluence in $\text{e}^-/\text{\AA}^2$

3 Beam block

A Beam block must begin with the keyword **BEAM**. At least the **TYPE** must be specified. Depending on the chosen **TYPE**, further declarations may be required.

3.1 TYPE

With the keyword **TYPE**, the underlying beam implementation is chosen. Currently three distinct beam implementations exist:

TYPE TOPHAT defines a beam with uniform flux.

TYPE GAUSSIAN defines a beam with a 2-dimensional Gaussian flux profile. The full-width half-maximum must be specified with the **FWHM** keyword (see section 3.3).

TYPE EXPERIMENTALPGM defines an arbitrary beam profile as defined in a given portable graymap file (.pgm) specified with the **FILE** keyword (see section 3.6). The pixel size of the image file also has to be specified with the **PIXELSIZE** keyword (see section 3.7).

3.2 FLUX

FLUX F specifies the total beam flux in photons per second. The flux parameter F can be specified in scientific notation (e.g. **1.3e12**). Note: The flux that should be specified is the flux incident on the crystal after collimation.

3.3 FWHM

FWHM X Y

This defines the X and Y FWHM of the beam respectively in the RADDOSE-3D coordinate system. For a horizontal goniometer, X = vertical and Y = horizontal, with it being the opposite for a vertical goniometer. FWHM is not needed if a Top-Hat beam is used.

3.4 ENERGY

ENERGY F

ENERGY F KEV

For a monochromatic beam, the energy specifies the incident photon energy in keV. The optional keyword **KEV** can be appended for human readability of the input file.

If the beam is polychromatic (as indicated by the specification of an **ENERGYFWHM**, see section 3.5) the energy specifies the mean of the distribution of photon energies, which is assumed to be a Gaussian distribution.

3.5 ENERGYFWHM

ENERGYFWHM F

specifies the FWHM of the distribution of incident photon energies in keV, which is assumed to be a Gaussian distribution. Not needed if the beam is monochromatic.

3.6 FILE

FILE *SOMEPGMFILE.PGM*

specifies the location of the .pgm file that describes the beam profile. NOTE: the .pgm file should be in ASCII format. The pixel size of the image file also has to be specified with the **PIXELSIZE** keyword (see section 3.7).

3.7 PIXELSIZE

PIXELSIZE *X Y*.

X and *Y* specify the horizontal and vertical size of the pixels in micrometers in the .pgm file respectively.

3.8 COLLIMATION

COLLIMATION RECTANGULAR *X Y*

specifies the *X* and *Y* collimation of the beam respectively in the RADDOS-3D coordinate system. For a horizontal goniometer, *X* = vertical and *Y* = horizontal, with it being the opposite for a vertical goniometer. Collimation delimits where the beam has non-zero intensity.

COLLIMATION CIRCULAR *X Y*

specifies the *X* and *Y* collimation of a circular or elliptical beam respectively in the RADDOS-3D coordinate system. *X* and *Y* are the diameter of the circle/ellipse defined by the slits. For a horizontal goniometer, *X* = vertical and *Y* = horizontal, with it being the opposite for a vertical goniometer.

3.9 PULSEENERGY

PULSEENERGY *F*

specifies the energy of the XFEL pulse in mJ. To be used instead of FLUX when the SUBPROGRAM XFEL input is used.

4 Wedge block

A Wedge block must begin with the keyword **WEDGE**.

WEDGE *A B*

A and *B* define the start and end angle of the rotation in degrees ($^{\circ}$). At 0° the front face of the crystal (*X*-*Y* plane) is normal to the beam. Rotation is right handed about the *Y* axis, as shown in figure 3).

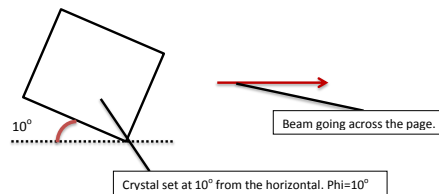


Figure 3: Schematic of angles for **WEDGE**. Figure courtesy of John Bremridge.

4.1 EXPOSURETIME

EXPOSURETIME *F*

specifies the total exposure time for this wedge in seconds.

4.2 ANGULARRESOLUTION

ANGULARRESOLUTION *F*

specifies the angular step size used for wedge iterations in degrees ($^{\circ}$). Defaults to 2° .

Note: If very small wedges are being used e.g. $< 5^{\circ}$ then the angular resolution should be decreased.

4.3 STARTOFFSET

STARTOFFSET *X Y Z*

offset translation in μm applied to the crystal relative to the origin (defined as the intersection of the beam and the aligned goniometer axis) for the starting position of the wedge. Defaults to 0 0 0.

4.4 TRANSLATEPERDEGREE

TRANSLATEPERDEGREE *X Y Z*

translation of the goniometer during exposure in $\mu\text{m}/^{\circ}$ for helical scanning, leading to improvements in dose distribution. Defaults to 0 0 0.

4.5 ROTAXBEAMOFFSET

ROTAXBEAMOFFSET *F*

the offset in μm along *X* (vertical in most set-ups) between the beam axis and the rotation axis. Used to create 'offset' scanning for improvements in dose distribution. Defaults to 0 μm .